TNT 919

Sample Blank Vials


for sample blank correction


Special Notes (For more detailed information: HACH Procedure Manual)

- Please read **Safety Advice** and **Expiration Date** on package.
- **Information:** TNT 846 (ortho-Phosphate) and TNT 848 (Aluminum) require the use of TNT 919, a special bar-coded Sample Blank Vial.
- **For special attention:** Please note that TNT 919 is not applicable for all other TNTplus™ vial tests (see Table 1).


Sample Blank Procedure for TNT 846 (ortho-Phosphate) and TNT 848 (Aluminum) with TNT 919

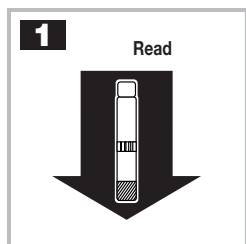Prepare a sample vial (TNT 846 or TNT 848) following the working procedure. Insert the sample vial (TNT 846 or TNT 848) into the photometer. This will launch the correct method and display the **uncorrected result**.

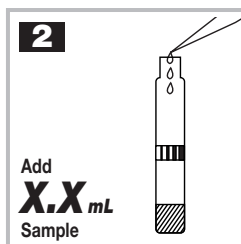Take a TNT 919 vial to prepare a Sample Blank for TNT 846 or TNT 848, considering the modifications according to Table 1. Insert the TNT 919 Sample Blank Vial with your prepared sample into the instrument. **The barcode tells the instrument that this is the Sample Blank.**

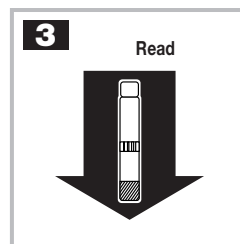If the Sample Blank value falls within the allowable range, this value will be used to correct the result automatically. The instrument subtracts the Sample Blank from the uncorrected result.

Summary	What is a Sample Blank
Colored and turbid samples can cause high-bias results to be obtained from photometric analyses. The influence of colors and/or turbidity can be eliminated by determining the sample blank.	A sample blank is obtained by performing a vial test without the reagent responsible for the color formation reaction . The result is attributable purely to the characteristics of the sample (color/turbidity). The sample blank can be deducted from the uncorrected result that is obtained when performing a vial test with a colored and/or turbid sample.

Table 1

Sample blank procedure							
Parameter	Test # TNTplus™	Vial	sample volume	dist Water volume	cap	reagent A	reagent B
Nutrients							
Ammonia	TNT 832	test Vial	0.2 mL		Zip (side w/out reagent)		
	TNT 830	test Vial	5.0 mL		Zip (side w/out reagent)		
	TNT 831	test Vial	0.5 mL		Zip (side w/out reagent)		
Nitrate	TNT 835	test Vial	1.0 mL	0.2 mL	original cap		
	TNT 836	test Vial	0.2 mL	1.0 mL	original cap		
Nitrite	TNT 839	test Vial	2.0 mL		Zip (side w/out reagent)		
ortho-Phosphate	TNT 844	test Vial	0.5 mL		Zip (side w/out reagent)		0.2 mL
	TNT 843	test Vial	2.0 mL		Zip (side w/out reagent)		0.2 mL
	TNT 845	test Vial	0.4 mL		Zip (side w/out reagent)		0.5 mL
	TNT 846	TNT 919	5.0 mL	1.0 mL	red stopper		
Metals (without digestion)							
Aluminum	TNT 848	TNT 919	3.0 mL		red stopper	2.0 mL	1 spoon full
Chromium (VI)	TNT 854	test Vial	2.0 mL		Zip (side w/out reagent)		
Nickel	TNT 856	test Vial	2.0 mL	0.2 mL	red stopper		

Note: For more detailed information see the HACH Procedure Manual.

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In the U.S.A. – Call toll-free 800-227-4224
Outside the U.S.A. – Contact the HACH office or distributor serving you.